


EVOLUTION[®]

PATIENT POSITIONERS

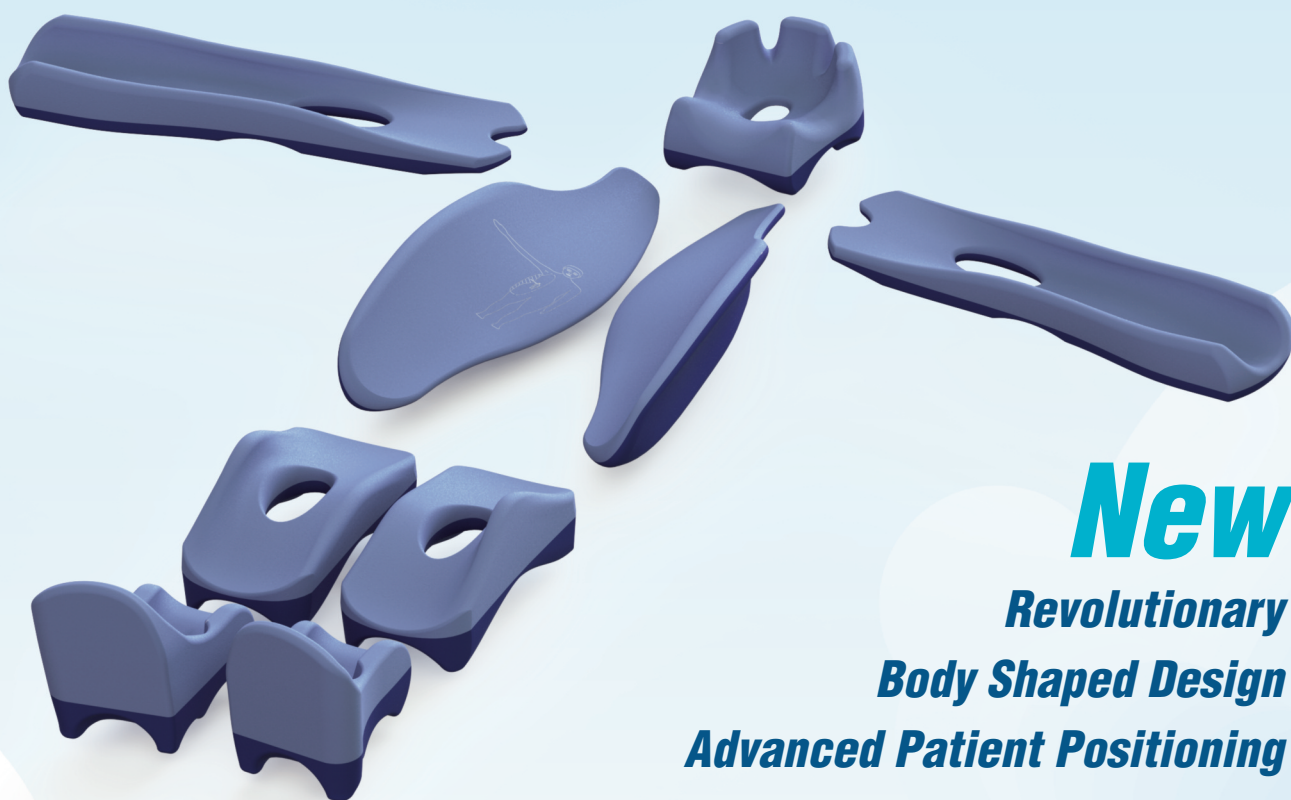

New
Revolutionary
Body Shaped Design
Advanced Patient Positioning
Relieves Pressure, Friction & Moisture
Cuts Risk of Decubitus Ulcers
Burst & Tear Resistant
Easily Cleaned

Just as the gel and foam devices of yesteryear represented an important evolutionary transition away from sand filled positioners, the EVOLUTION range of patient positioning devices represents the next important technological leap forward.

View our Evolution Product Video at
www.vascocare.com


VASCO CARE
MEDICAL LIMITED



Clinical Tests prove that the EVOLUTION Head Positioner provides superior pressure relief compared to what is currently on the market. VascoCare is presenting the EVOLUTION Supine Head Positioner, which was developed at our research laboratory using our new patented technology that is unique to the market.

Enhanced Head Stability and Patient Safety

Allows ease of Endotracheal Intubation



Style 1



Style 2



Redistributes pressure over a greater area of the patients head and neck

Helps prevent the formation of Decubitus Ulcers


265lb Person on Evolution Patient Positioner Large Head Style 1 SHP14

0 5 15 25 35 45 55 65 75 85 95 mmHg

Zero Pressure on the Occipital area helps to prevent Decubitus Ulcers, Pressure Alopecia and Fluid Pooling

Fingered Lateral Supports provides excellent pressure redistribution and stable head support with comfortable fit

Raised neck area supporting cervical vertebrae and trapezoid muscle

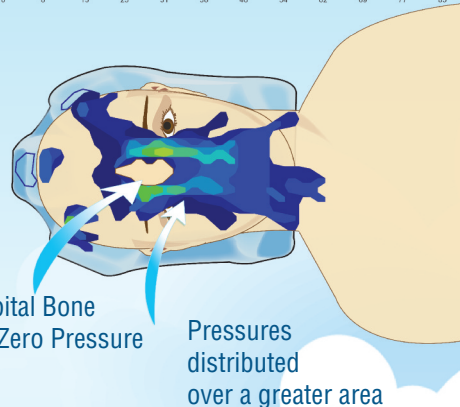

Occipital Bone with Zero Pressure

Pressures distributed over a greater area

Smooth contoured surfaces allows maximum patient immersion and ease of cleaning

Air Ventilation through the base maintains normothermia at the occipital bone

Recessed handles for patient manipulation, helps eliminate risk of injury

Product Information permanently engraved into coating for regulatory and traceability requirements

Concave base aids self adjusting motion and allows air ventilation

Non-Slip base grips to operating table mattress

Covers:

REF CSHP12

REF CSHP14

REF CSHP23

REF CSHP24

Supine Head Positioners:

SHP12 Small

SHP14 Large

SHP23 Medium

SHP24 Large



Style 1



Style 1




Style 2



Style 2



Enterprise and Research Incubation Centre
Institute of Technology Carlow, Kilkenny Road, Carlow, Ireland
T: 00353 (0)59 9175260 E: customer@vascocare.com www.vascocare.com
© Copyright VascoCare Medical Ltd 2/2012
REF: Flyer/SHP/004

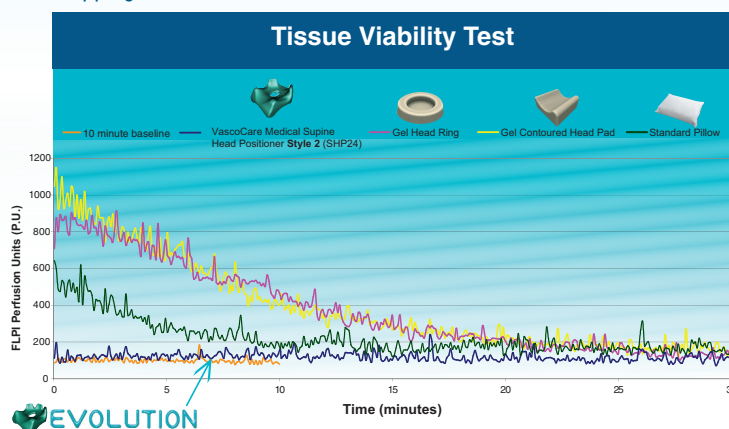

Authorised Distributor

